

Autoclave Tape



Cost Effective

- Conforms to ISO 11140-1 Class 1
- Easy to interpret
- Cost effective
- Individual rolls marked with lot number and expiry date
- Non-toxic
- Lead free


Autoclave Tape

Autoclave tape for steam sterilization enables the user to see at a glance that an item has been processed. The easy to interpret black stripes prove that the pack has been exposed to the steam sterilization process.

The tape securely sticks to non woven, paper, board, metal, plastic and glass and is manufactured from treated crepe paper and coated with high performance cross linked heat sensitive adhesive

Order Details

V3504040	Autoclave Tape	19mm 48 rolls per bag
V350 4045	Autoclave Tape	24mm 36 rolls per bag

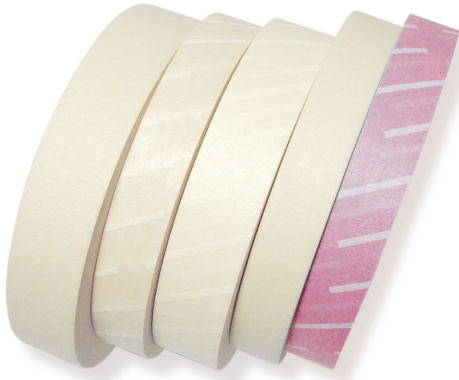

V3504040 19mm x 50m roll

V3504045 24mm x 50m roll

Name:	Autoclave tape for steam sterilization
Classification:	EN ISO 11140-1 Class 1
Operating temp:	121°C – 134°C
Time	As per specific sterilization cycles minimum 134C 3.5 minutes or equivalent
Storage conditions:	10-30°C
RH:	30-70%
Lot no:	yes on core
Prod date:	yes on outer carton
Exp date:	yes on core
Correct symbols:	yes in accordance with the standard
Indicator:	yes
Expiry date:	2 years
Extra Information:	
Substrate	Medical crepe paper
coated	yes
Recordable Information:	Information can be written on the tape with a Valisafe PIP Colour change: white to black
Packaging:	24mm width tape = bags of 36 rolls 19mm width tape = bags of 48 rolls